# Zen II

PLUS

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## *DearMiss* DWS (DEEP WRINKLE SOLUTION)

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## *DearMiss* SWS (SHINING WHITE SOLUTION)

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## *DearMiss* PCS (PROBLEM CONTROL SOLUTION)

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# TABLE OF CONTENTS

INTRODUCTION .....	P24
SAFETY PRECAUTIONS .....	P24~25
PRECAUTIONS BEFORE USE .....	P26
PRECAUTIONS IN USE .....	P27
PRECAUTIONS OF MAINTENANCE AND INSPECTION .....	P28
NAMES OF COMPONENTS .....	P29~31
HOW TO USE .....	P32~34
HOW TO APPLY BEAUTY ESSENCE AND GEL .....	P34
PREPARATION BEFORE TREATMENT .....	P34
TREATMENT PROCEDURE .....	P35~38
Dr.Arrivo Zeus IV DEDICATED BEAUTY ESSENCE .....	P39~41
PRODUCT SPECIFICATIONS .....	P42

# INTRODUCTION

- ◆ Thank you for purchasing ARTISTIC&CO. product.  
In order to use the device properly and safely, it's necessary to read this user manual carefully and follow all safety precautions. Please include this user manual if you are borrowing or giving to another person.
- ◆ This instruction manual contains the warranty card.  
Confirm "Date of purchase" and "Purchase shop", etc. are entered on the warranty and store it carefully.

## SAFETY PRECAUTIONS

Be sure to read this user manual carefully before using, always observe the following safety precautions.

To reduce the risk of injury, loss of life, and damage to any property, the followings symbols are used to classify and describe the type of instructions to be observed. The marks and meanings of symbols are shown as below.

	<b>DANGER</b>	Failure to observe DANGER signs will result in hazard to personnel is high likely to result in serious injury or death.
	<b>WARNING</b>	Failure to observe WARNING signs could result in hazard to personnel possibly resulting in serious injury or death.
	<b>CAUTION</b>	Failure to observe CAUTION signs could result in injury to personnel and damage to the device and property. ※ Property damage includes damage to house, household belongings, domestic animals and pets.

### ■ MEANINGS OF SYMBOLS

	This symbol indicates prohibition. ※ inside or near it is indicated a special instruction.
	This symbol indicates requirement that must be followed.
	This symbol indicates caution.

# SAFETY PRECAUTIONS

## WARNING

### To avoid accidents and physically ill.

- ◆ Do not use this device with any of the following medical electrical equipment. Malfunctions of the medical electrical equipment may cause accidents or illness.

Prohibition

- Implantable medical electrical equipment such as pacemaker.
- Implantable medical electrical equipment such as artificial heart-lung device
- Wearable medical electrical equipment such as electrocardiograph or dropper.

Prohibition

- The following persons should not use the device.
  - Person who has heart disease, abnormal blood pressure or malignant tumors.
  - Person who has disease with a fever, infectious diseases or acute diseases.
  - Person who is undergoing dental treatment.
  - Person who has weak immune systems or uses immunosuppressants due to immunity disorders like AIDS etc.
  - Person who has hemophilia, bleeding disorders, diabetes or facial neuralgia.
  - Person who has paresthesia or dysalgia.
  - Person who has keloid diathesis, facial melanosis, photosensitivity or pigmentation disorders.
  - Person who has capillary dilatation due to long-term steroidal hormone treatments or liver dysfunction.
  - Person who has allergies, atopic dermatitis, skin diseases, sensitive skin or weak skin.
  - Person who has medical history of epilepsy or gets sick due to flashing lights.
  - Person who has sequela of head or ophthalmologic sequela due to traffic accidents etc.
  - Person who had eye operations (in the case of LASIK, consult with the doctor before using).
  - Person who is pregnant, on period or breast-feeding
  - Children or persons who can not express themselves.
  - Person who is in a drunken or hangover state.
  - Person who is taking sleeping medicine.
  - Person who is under state of extreme fatigue or physical discomfort.
  - Person who had skin problems such as spots or rashed due to artificial fibers.
  - Person who has or had skin inflammation before due to cosmetics.
  - Person who recently had chemical peeling treatment at salons etc.
  - Person who had skin problems after using beauty device before.

Prohibition

- Do not use the device on the following areas.
  - Except specified parts.
  - Around heart, eyeballs, upper eyelids, lips and musous parts such as inside the mouth.
  - Pathological spots, moles or viral warts.
  - Areas where skin is irritated or inflamed due to excessive sunburn.
  - Areas where skin suffers from Cuts, eczema, psoriasis, herpes, puerlent acnes or pimples.
  - Areas where skin suffers from pathological spots due to itching, flushing, physical stimulus etc.
  - Cuts where surgery was performed, parts, anus, areas around adam's apple or ears.
  - Any areas where has performed hair removal or shaving recently.

※ If you are seeing a doctor, consult with the doctor before use.

# PRECAUTIONS BEFORE USE

## ⚠ WARNING

 <p>Prohibition</p>	<ul style="list-style-type: none"> <li>● Do not use this device with other beauty device at same time.</li> <li>● Do not use this device for any purpose other than beauty and under conditions of malfunction or damage. This device is a equipment specialized for use in beauty. Do not use for therapeutic or any other purposes.</li> <li>● Do not stare or look directly at the LED light.</li> <li>● Do not use the device continuously on the same areas. Failure to comply may cause or aggravate skin damage, infection, inflammation, bleeding, or pigmentation.</li> <li>● Do not use on pets.</li> </ul>
 <p>Prohibition of use in a wet area</p>	<ul style="list-style-type: none"> <li>● Do not use or store the device in the bathroom or anywhere with high humidity. Do not leave the device where is easily splashed. (such as bathtub, washbashi.)</li> </ul>
 <p>Prohibition of getting wet</p>	<ul style="list-style-type: none"> <li>● Do not immerse the device in water or pour or splash water on it. (Failure to comply may cause electric shock fire or malfunction.)</li> </ul>
 <p>Prohibition of wet hands</p>	<ul style="list-style-type: none"> <li>● Do not touch the device with wet hands. (Failure to comply may cause electric shock or injury.)</li> </ul>
 <p>Prohibition</p>	<ul style="list-style-type: none"> <li>● Do not use the device near a flammable substances such as thinner, benzene, polish remover or spray.</li> <li>● Do not use with anything other than the dedicated adapter. Do not connect the dedicated adapter to other equipments.</li> </ul>
 <p>Unplug</p>	<ul style="list-style-type: none"> <li>● Be sure to hold the front of the plug when plugging in or out the adapter. Do not grab the cord to plug in or out the adapter.</li> <li>● Do not use the adapter when it is damaged or loose.</li> </ul>
 <p>Prohibition</p>	<ul style="list-style-type: none"> <li>● Do not damage the adapter cord like wrecking, twisting, pulling, placing heavy items on or squeezing it.</li> <li>● Do not overload the outlet. (Failure to comply may cause electric heat or fire.)</li> </ul>
 <p>Compulsion</p>	<ul style="list-style-type: none"> <li>● Make sure to use household outlets.</li> <li>● Make sure to insert plug firmly to adapter.</li> </ul>

# PRECAUTIONS IN USE

## ⚠ WARNING

 Unplug	<ul style="list-style-type: none"> <li>● Stop using the device, turn off the power and disconnect the AC adapter immediately in the following conditions.             <ul style="list-style-type: none"> <li>• Earthquakes or lightning.</li> <li>• Smoke, excess heat, odor or distortion from the device, adapter or cord.</li> </ul> </li> </ul>
 Prohibition	<ul style="list-style-type: none"> <li>● During treatment, this device might get warm but there will not be a safety issue.</li> <li>● Do not contact the device with any metal or electric appliances.             <ul style="list-style-type: none"> <li>※ mobile phone, television, etc.</li> </ul> </li> </ul>
 Prohibition of disassembly	<ul style="list-style-type: none"> <li>● Do not disassemble, modify, or try to repair the device yourself. For repair, please inquire at dealer.</li> </ul>
 Compulsion	<ul style="list-style-type: none"> <li>● Make sure to remove all the metal accessories, glasses and contact lenses during treatment.             <ul style="list-style-type: none"> <li>※ metal accessories: pierced earrings, necklaces, watches, etc.</li> </ul> </li> <li>● Make sure to use only level 1 during your first treatment and until you are accustomed to use.</li> <li>● Do not press the device head to the skin strongly or use it on the same spot over 3 seconds. Failure to comply may cause skin problems.</li> <li>● Keep eyes closed during treatment around eyes.</li> <li>● If you feel pain on skin or body, or any abnormalities, discontinue use immediately and consult a doctor.</li> </ul>

## ⚠ CAUTION

Precautions of use of main body      Read the following sections carefully before use.

1	<p>Perform a special treatment to each part once per week.</p> <ul style="list-style-type: none"> <li>● Perform a special treatment twice a week in the first and second weeks. From the third week, perform a treatment about once a week according to your skin conditions.</li> </ul> <p>Perform a daily treatment for each part every day or every two days.</p> <ul style="list-style-type: none"> <li>● Use the device according to your skin conditions.</li> </ul>
2	<ul style="list-style-type: none"> <li>● Before using the device for a facial treatment, remove makeup and wash face to clean up impurities such as makeup cosmetics and sebum, and apply the dedicated beauty essence or gel.</li> </ul>
3	<p>Skin conditions are different depending on physical conditions.</p> <ul style="list-style-type: none"> <li>● During treatment, carefully observe your skin conditions, if your skin shows any abnormalities like rashes, redness or itchiness, discontinue use immediately and consult a doctor.</li> </ul>

# PRECAUTIONS OF MAINTENANCE AND INSPECTION

## ⚠ WARNING

 Unplug	<ul style="list-style-type: none"> <li>● Unplug the adapter from the outlet when cleaning the main body. (May cause fire or an accident.)</li> </ul>
 Prohibition	<ul style="list-style-type: none"> <li>● Do not allow adapter or plug to contact with pins or dust. (Failure to comply may cause fire, electric shock or damage.)</li> <li>● Do not drop or hit the device.</li> <li>● Do not wrap the code around the adapter when storing the device</li> </ul>
 Unplug	<ul style="list-style-type: none"> <li>● Make sure to turn off the power, and disconnect the AC adapter with the outlet when not using. (Failure to comply may cause fire or damage.)</li> </ul>
 Prohibition	<ul style="list-style-type: none"> <li>● Keep away from children and do not allow children to play with the device or attachments.</li> <li>● Do not leave the device in the following locations. <ul style="list-style-type: none"> <li>• A location where is exposed to direct sunlight, or high humidity, or dust.</li> <li>• A location with high temperature or subzero temperature, such like cars under the sun.</li> <li>• A location where is flammable, or outdoor with wind and rain.</li> <li>• A location where there is chemicals, gas, strong vibration.</li> <li>• A location where there is near a heater, heating device, or magnet</li> </ul> </li> </ul>
 Compulsion	<ul style="list-style-type: none"> <li>● How to clean <ul style="list-style-type: none"> <li>• Wipe the device by well-wrung soft cloth. (If using hard cloth may cause damage to the device, or deterioration in functions)</li> </ul> </li> <li>※ Do not wipe the device with solvents such as alkaline agents, alcohol, thinner or benzine, polish remover, polishing power, hot water or chemicals.</li> <li>※ Do not leave the device with the essence or gel on it after use.</li> </ul>

## Inspection

- If you have not use the device for longer than 6 month, make sure to check there are any abnormalities (for example, the state of LED lighting, or whether the cracks), before starting use again. Use the device under malfunction may cause an accident.

# NAMES OF COMPONENTS

## PRODUCT STRUCTURE

## Accessories

- AC adapter
- Stand
- User Registration Card
- Guarantee Card

# NAMES OF COMPONENTS

## 3D movable head

A unique device head mechanism can flex in all directions for movement.

Follow the face curves, close to skin and perfect fitted.


Stress-free operation is available.

## 24KEp electrode

The surface is surrounded with 24-karat gold which is more gentle on the skin.

※ This metal is really soft and easily scratched. Please handle it with sufficient care. If there are any small scratches, it will not any problem with quality of the device. It is safe to use.
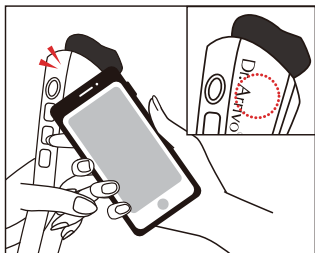
## NFC tag-equipped model

Available to check the serial number of each device by placing the NFC-compatible device over this device.

※ Follow the figure left.

### Compatible models





iOS: iPhone XS/XR or newer  
Android: all

※ Models iPhone 5S/5C to iPhone X require a compatible app.

(Search for "NFC Reader" in the App store)

※ Not available with NFC-incompatible models.

# HOW TO USE

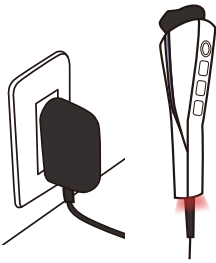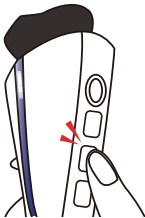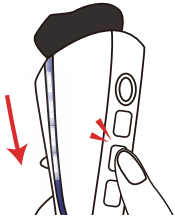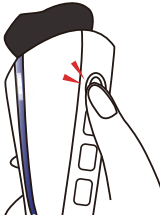
1	<p><b>Plug the AC adapter into a outlet and connect to the main device.</b></p> <ul style="list-style-type: none"> <li>◆ Each lamp is lit with a starting sound.</li> <li>※ When using this product for the first time, remove the protective sticker (the transparent sticker) from the elements.</li> </ul>	
2	<p><b>Select the mode.</b></p> <ul style="list-style-type: none"> <li>◆ Select from the following 3 modes: RF mode, M/U mode, and EMS mode.</li> <li>※ &lt;RF&gt; button: RF mode</li> <li>※ &lt;M/U&gt; button: M/U mode</li> <li>※ &lt;EMS&gt; button: EMS mode</li> <li>◆ The color of the side lamp changes depending on the selected mode.</li> <li>※ To change the mode, press the button of the mode being used to turn it off and then select the new mode.</li> </ul>	
3	<p><b>Select the level.</b></p> <ul style="list-style-type: none"> <li>◆ Press the button of the selected mode to select the level.</li> <li>※ It's recommended to use output level 1 setting during your first treatment and then adjust to the output level which is comfortable for you.</li> <li>※ To reduce the level, press the switch of the mode being used to turn it off and then reset the new level again.</li> <li>※ See P. 31 for more information.</li> </ul>	
4	<p><b>Select the LED.</b></p> <ul style="list-style-type: none"> <li>◆ Press the &lt;LED&gt; button when using the RF mode, the M/U mode, or the EMS mode with 5 types of LED modes can be selected.</li> <li>※ See P. 32 for more information.</li> </ul>	

# HOW TO USE

## Name of parts

## Description of level lamp
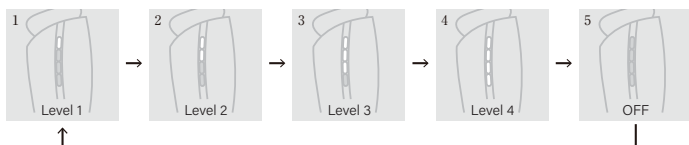
## Functional description

Button	Mode	Function
RF	RF mode	Warm skin gently from the deep to lead to skin with resilience.
M/U	M/U mode	MFIP mode and UP mode are out putting alternately. MFIP(medium frequency intermittent pulse) =ULTRA PULSE outputs intermittently. ULTRA PULSE=Insert pulse/Mid-high frequency/EMS outputs simultaneously.
EMS	EMS mode	Output a pulse specialized in working muscles.

## Description of machine sound

During use, a <Beep> sound can be heard.	When turning on a mode and use it, you can hear the sound around once every 10 seconds.
	The sound can be heard once every second from before around 5 seconds of a stop. (The power automatically off after 10 minutes have passed from the start of a treatment.)


# HOW TO USE

## Description of LED

- ※Do not look directly at the LED light.
- ※During the treatment around eyes, do not use the device with eyes open.
- ※The treatment can be performed while the LED light is off too.

<div>Press once All colors of the LEDs</div> <div></div>	3 different types of LED lights are installed, red, green, and blue. Recommended for who has several skin troubles.
<div>Press twice Blue LED 415nm(±10nm)</div> <div></div>	Lead to healthy skin.
<div>Press three times Green LED 530nm(±10nm)</div> <div></div>	Lead to translucent and impressive skin.
<div>Press four times Red LED 620nm(±10nm)</div> <div></div>	Lead to elastic and lustrous skin.
<div>Press five times Flashing LED</div> <div></div>	Recommended to use in the final process of the treatment.

# HOW TO APPLY BEAUTY ESSENCE AND GEL

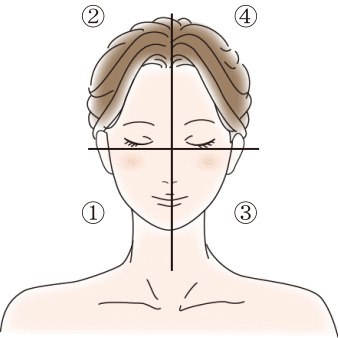When performing a treatment for a full face, divide your face into four areas and apply beauty assence or gel on each treatment area in turn.

## PREPARATION BEFORE TREATMENT

- ① Remove accessories such as necklaces or rings and all metal objects that you are wearing.
- ② Remove makeup and wash your face to remove impurities such as oil and dirt from your skin.
- ③ According to the mode being used, apply the beauty essence or gel on the treatment area.

Dedicated beauty essence	Dedicated gel
<div></div> <div>PE GOLDEN BEAUTY THE SERUM</div> <div>Mode in use</div> <div>M/U mode</div>	<div></div> <div>PE GOLDEN BEAUTY THE MASSAGE GEL</div> <div>Mode in use</div> <div>RF mode EMS mode</div>

# TREATMENT PROCEDURE

## Special treatment once a week

◆ RF mode : 10 minutes + M/U mode : 10 minutes = Total 20 minutes

- ① Apply the gel on the treatment area when using the device in the RF mode.
  - ② Apply the beauty essence on the treatment area when using the device in M/U mode.
- ※ Do not put the device on the same area of your skin for longer than 3 seconds to prevent skin problems.
- ※ Treatment only either RF mode or M/U mode can be performed for 10 minutes.

	<p>Move the device from under the chin to the below ear and then down to the shoulder.</p>
	<p>Move the device from under the mouth to the temple. Next, move the device from the nostril to the temple. Finally, move the device from under your eye to the temple.</p>
	<p>Move the device from the eyebrows to the hairline.</p>

- ★ Perform each procedure from ① to ③ in the illustration two or three times.
- ★ Move the device constantly and slowly and do not put the device on same spot for a long time.
- ★ When the beauty essence or gel becomes dry, add more.

# TREATMENT PROCEDURE

## Daily treatment

### ◆EMS mode = 10 minutes

① Apply the gel on the treatment area and then perform a treatment.

※ During the treatment, apply the gel more if the treatment area becomes dry.


When the treatment area becomes dry, it may cause irritation on your skin and the reduction of its effect.

※ All nine electrode of the surface should contact the treatment area closely. (Do not make space between the all nine electrode of the surface and the treatment area.)

※ Do not use the device many times a day as it may hurt your skin.

※ During the treatment, if the surface becomes rough, please check the device and clean up with soft cloth or tissues.

※ Do not use the device if you are physically unwell.

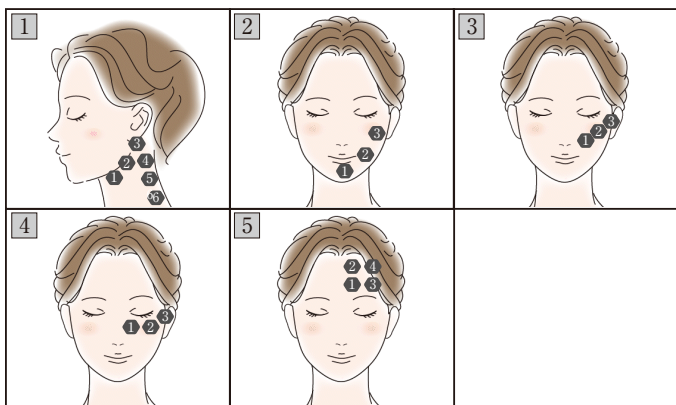★ Perform each procedure from 1 to 5 in the illustration two or three times.

★ Use the device on each treatment point for approximately 3 seconds and then move to the next point.

※ The treatment procedure in the illustration is a reference.

★ When the treatment area becomes dry, apply the gel or beauty essence more.

## WARNING

Do not use on the mucosa, the upper eyelids, the eyeballs, and around the ears and the Adam's apple

# TREATMENT PROCEDURE

## Body treatment


※ Perform a treatment on the clean skin.

※ Apply the gel or beauty essence on the treatment area.

### Upper arm

■ Special treatment : each arm 20 minutes

■ Daily treatment : each arm 10 minutes

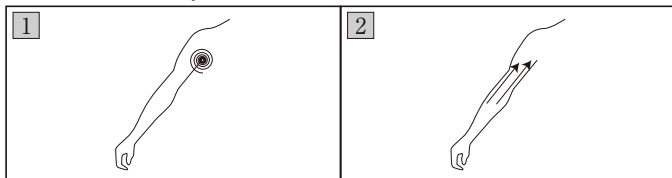
1 Move the device in small circles around the shoulder joint 2 or 3 times.

2 Move the device straight from around the elbow to the shoulder joint.

### Waist

■ Special treatment : 20 minutes

■ Daily treatment : 10 minutes

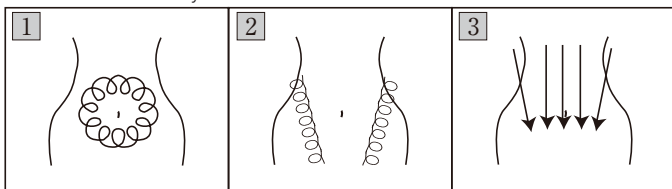1 Move the device in small circles around the navel.


2 Move the device in small circles from the side of the waist to the inguinal lymph node.

3 Move the device straight from the top to bottom.

### Hip

■ Special treatment : 20 minutes

■ Daily treatment : 10 minutes

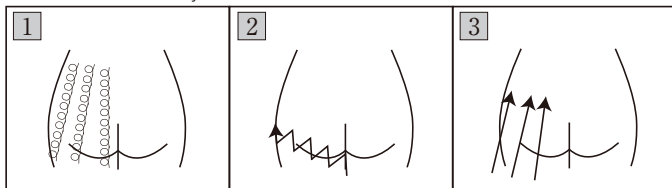1 Move the device in small circles from the groin to the lower back.

2 Move the device in zig-zags around the groin.


3 Move the device straight from the groin to the lower back.

# TREATMENT PROCEDURE

## Thigh (front and back)

■ Special treatment : each side of thigh 20 minutes

■ Daily treatment : each side of thigh 10 minutes

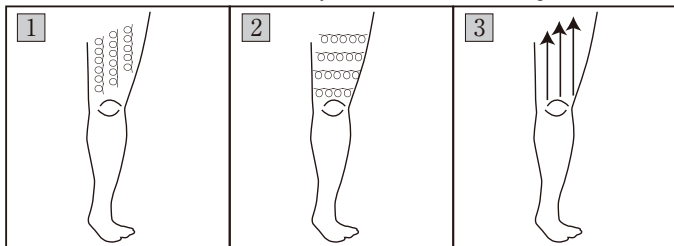1 Move the device in small circles from above the knee (the back of the knee) to the groin.

2 Move the device in small circles from side to side from above the knee (the back of the knee) to the groin.


3 Move the device straight from above the knee (the back of the knee) to the groin.

※ Perform in the same way for the front and the back of the legs.

## Calf

■ Special treatment : each calf 20 minutes

■ Daily treatment : each calf 10 minutes

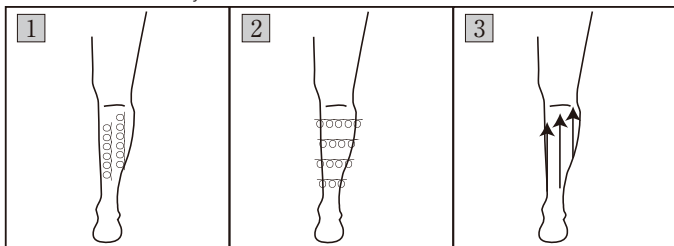1 Move the device in small circles from around the ankle to the back of the knee.

2 Move the device in small circles from side to side from above the ankle to the back of the knee.

3 Move the device straight from above the ankle to the back of the knee.

## After treatment

① After the treatment, gently massage any remaining serum into the skin.

If there is any residual gel and you are concerned about stickiness,

② After the treatment, wipe away the device with a soft cloth or tissues. After that, pat the device dry with tissues.

# Dr.Arrivo Zeus II PLUS

## DEDICA BEAUTY ESSENCE

### PE GOLDEN BEAUTY THE SERUM

#### Beauty essence for " Arrivo series " facial treatment

A dedicated beauty essence pursuing the affinity between Ultra Pulse mode and MFIP mode. Leads to resilient skin that is lifted.

Approaches to moisture and resilience required for maintaining beautiful and healthy skin. Conditions mature skin. Additionally, richly formulated with the moisturizing ingredients, achieving the texture that has a good affinity with the dedicated beauty equipment, " Arrivo series ".

Use in the M/U mode.

PE Golden Beauty the Serum    40mL JPY15,400 (tax included)  
120mL JPY39,600 (tax included)

### PE GOLDEN BEAUTY THE BODY

#### Beauty essence dedicated for " Arrivo series " body treatment

Formulated with originally-blended beauty ingredients for toning the skin.

Additionally, you can enjoy a pleasant treatment time with the scent of carefully selected herbal extracts. A massage gel which also focuses on a texture that spreads well and blends into the skin quickly.

Use in the RF mode, the M/U mode or the EMS mode.

PE Golden Beauty the Body    200g JPY16,500 (tax included)

# Dr.Arrivo Zeus II PLUS

## DEDICA BEAUTY ESSENCE

### PE GOLDEN BEAUTY THE MASSAGE GEL

#### Massage gel dedicated for " Arrivo series "

The moist veil adheres to the skin.  
Bring out brightness and transparency.

Water-soluble massage gel richly formulated with beauty essence ingredients.  
Leads to translucent and soft skin with continuous use. An original formula that focuses on moderate texture and excellent spread to easily perform a treatment.

※ After massage, it is recommended to wipe away gel left on the skin.

Use in the RF mode or the EMS mode.

PE Golden Beauty the Massage Gell 200g JPY5,500 ( tax included)

### PE THE DAY

#### Beauty essence for " Arrivo series " iontophoresis

An elegant drop of " THE DAY " . . .For sophisticated skin every day.

For a wide range of age groups and various skin types. Gently wraps the skin with moist veil to protect it.

Use in the M/U mode.

PE The Day 30mL JPY7,700 (tax included)

# Dr.Arrivo Zeus II PLUS

## DEDICA BEAUTY ESSENCE

### *Dear Miss* DWS (DEEP WRINKLE SOLUTION)

Intensive beauty essence " DWS " dedicated for red LED treatment.

Lead to lustrous and elastic skin.

---

Dear Miss DWS    30mL   JPY 16,500( tax included)

### *Dear Miss* SWS (SHINING WHITE SOLUTION)

Intensive beauty essence " SWS " dedicated for green LED treatment.

Lead to translucent nad impressive skin.

---

Dear Miss SWS    30mL   JPY 13,200( tax included)

### *Dear Miss* PCS (PROBLEM CONTROL SOLUTION)

Intensive beauty essence " PCS " dedicated for blue LED treatment.

Lead to healthy skin.

---

Dear Miss PCS    30mL   JPY 13,200( tax included)

# PRODUCT SPECIFICATIONS

	Product Name	Dr.Arrivo Zeus II PLUS
	Power Rating	Input AC 100-240V~50/60Hz    Output DC 12V • 1.5A
	Power Consumption	Maximum 17W
	Output Frequency	90kHz / 200kHz / 100Hz
	LED Wavelength	Blue415nm(±10nm)/Green530nm(±10nm)/Red620nm(±10nm)
	Size	Maximum 189mm×55mm×58mm
	Weight	approximately 151g
	Timer	10 minites

*Top of the Beauty*